



The future of vision

# OCULFIT CONICAL ORBITAL IMPLANT

Minor porous  
superfície




ANTERIOR ZONE



SIDE VIEW

# OCULFIT CONICAL Orbital Implant

|                             |  |           |           |
|-----------------------------|--|-----------|-----------|
| Brand name                  | OCULFIT  |           |           |
| References                  | 80016-POI  | 80018-POI | 80020-POI |
| Sphere diameter             | 16 mm  | 18 mm     | 20 mm     |
| Description and Indications | <p>OCULFIT CONICAL orbital implants are implants designed for orbital reconstruction following an evisceration or enucleation procedure.</p> <p>They are made from high-density polyethylene, which has been proven to be well tolerated by tissues and which provides stability, giving surgeons biocompatible porous materials for orbital reconstruction.</p> <p>The advantage of the structure of open, interconnected pores is that it enables vascularisation and the growth of soft tissue inside it. The anterior surface of Oculfit Orbital Implants is <b>smooth and porous with four shallow suturing tunnels in an aim to facilitate the attachment of the rectus muscles</b>. The smoother porous surface can help to minimise the exposure of the implants in the long term.</p> |           |           |
| Material                    | High density porous polyethylene.  |           |           |
| Sterilization method        | Ethylene Oxide (ETO).  |           |           |
| Supply                      | Sterile, single container. The devices are accompanied by an insertion accessory that facilitates their implantation, avoiding the adhesion of the tissues in its path.  |           |           |
| Shelf-life                  | 5 years  |           |           |

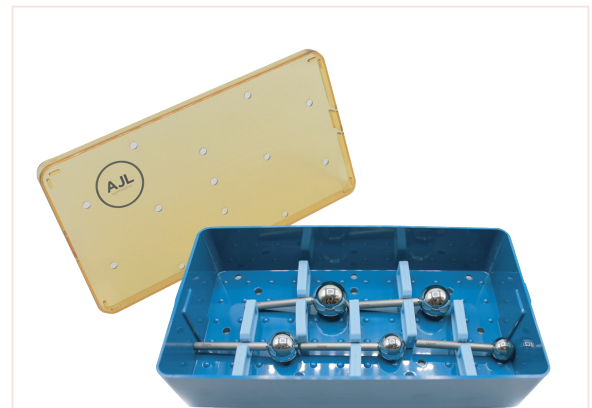

## PERFORATED CONFORMER



THE PERFORATED CONFORMER IS RECOMMENDED TO BE PLACED IMMEDIATELY AFTER THE SURGICAL PROCEDURE.

|            |                            |                             |                            |
|------------|----------------------------|-----------------------------|----------------------------|
| Brand name | AJL PERFORATED CONFORMER   |                             |                            |
| References | CF-S-AJL                   | CF-M-AJL                    | CF-L-AJL                   |
| Size       | Small:<br>22 x 19 x 1.2 mm | Medium:<br>24 x 21 x 1.2 mm | Large:<br>25 x 22 x 1.2 mm |

|            |                                    |
|------------|------------------------------------|
| Brand name | ORBITAL VOLUME SIZER SET           |
| Reference  | AMB-A1-SET                         |
| Size       | 14 mm, 16 mm, 18 mm, 20 mm, 22 mm. |

\* Product with CE marking

REV: 01. 03/2021